

Powered Instruments

for large bone surgery



de Souther
MEDICAL

powered instruments for large bone surgery

CONTENTS

MultiDrive™ system

Modular handpiece attachments	3
-------------------------------------	---

Electric powered

MDX MultiDrive™ battery and corded handpieces	5
---	---

Procedure specific battery and corded handpieces	6
--	---

Battery options	7
-----------------------	---

Battery charger and corded direct power	8
---	---

Air powered

MPX MultiDrive™ air powered handpieces	9
--	---

Procedure specific air powered handpieces	10
---	----

PLX Pulsed irrigation system	11
------------------------------------	----

Ancilliary equipment

Sterilisation containers	12
--------------------------------	----

Accessories	13
-------------------	----

All products in this catalogue are backed by the De Soutter Medical no nonsense warranty.
All instruments shown in this catalogue are CE marked.


The technologically advanced De Soutter Medical large bone surgical system sets new standards for orthopaedic power tools. Enjoy the versatility of real power, compactness and lightness with a comprehensive range of quick change attachments. Be it large bone or trauma, the De Soutter Medical large bone system covers them both with precision, performance and reliability.

The versatile large bone system



MODULAR SYSTEM

26 colour coded attachments with proven, quick change locking mechanism.
See page 3.

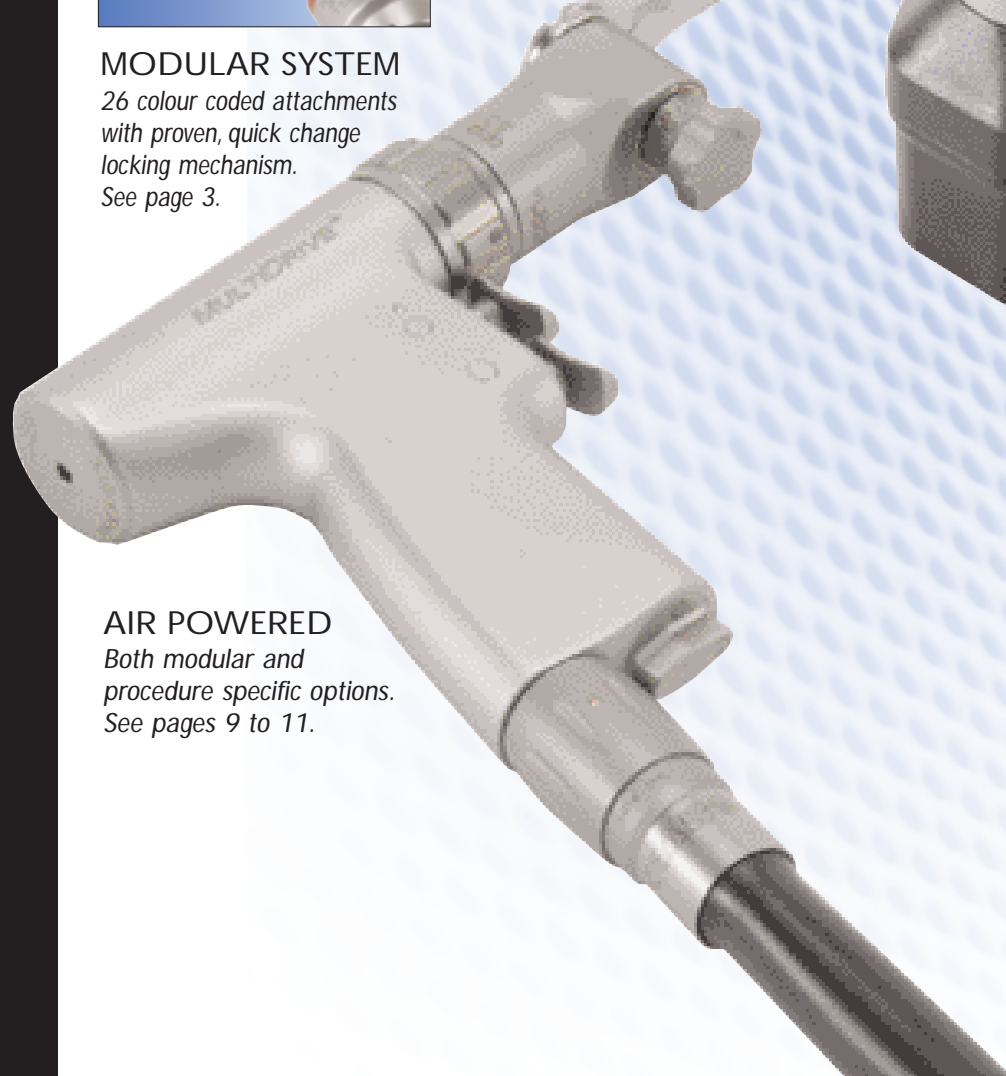

AIR POWERED

Both modular and procedure specific options.
See pages 9 to 11.

BATTERY AND CORDED

Both modular and procedure specific options with battery, aseptic battery and corded direct power.
See pages 5 to 8.



MultiDrive™ System

Modular handpiece and attachment system



MDX-610

Battery or corded
handpiece options.
See pages 5 - 8.



MPX-610

Air powered handpiece.
See page 9.

Drilling attachments



Order code	Chuck Description	Type	Speed rpm/cpm	Cannulation mm
11710	0 - 6.35 mm (1/4") Jacobs	DX-600	0 - 950	4.1
11720	0 - 4 mm (5/32") Jacobs	DX-600	0 - 950	3.1



11860	Hi-torque - 0 - 6.35 mm (1/4") Jacobs	DX-600	0 - 500	4.1
-------	--	--------	---------	-----



11740	0 - 3.2 mm (1/8") keyless	DX-600	0 - 950	2.4
11750	0.8 - 6.35 mm (1/4") keyless	DX-600	0 - 950	4.1

Quick release drilling attachments



11730	A.O. Synthes™/ASIF	DX-600	0 - 950	2.0
11770	Hudson	DX-600	0 - 950	4.1
11790	Trinkle	DX-600	0 - 950	4.1
12420	Zimmer	DX-600	0 - 950	4.1
14480	Hudson/Zimmer combi	DX-600	0 - 950	4.1



11760	Oscillating drill - 130° oscillation arc	OX-600	0 - 950	-
-------	---	--------	---------	---

See page 14 for range of drill bits
See Accessory catalogue No. 162 for range of taps.

Standard reaming attachments



11800	0 - 6.35 mm (1/4") Jacobs	RX-600	0 - 260	4.1
-------	---------------------------	--------	---------	-----



11810	Hudson	RX-600	0 - 260	4.1
11820	A.O. Synthes™/ASIF/PROTEK	RX-600	0 - 260	4.1
11830	Zimmer	RX-600	0 - 260	4.1
12900	Hudson/Zimmer combi	RX-600	0 - 260	4.1

Hi-torque reaming attachments



11850	Hudson/Zimmer combi	HX-600	0 - 150	-
11840	A.O. Synthes™/ASIF/PROTEK	HX-600	0 - 150	-

See page 13 for acetabular reamer shafts and shells


Sawing attachments




12390 ¹⁾	Reciprocating saw See Stericut™ catalogue for saw blade range.	CX-600	0 - 9000	
---------------------	--	--------	----------	--



12600 ²⁾	Reciprocating saw See Stericut™ catalogue for saw blade range.	CX-650	0 - 7500	
---------------------	--	--------	----------	--

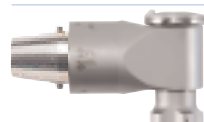

12380 ¹⁾	Large sagittal See Stericut™ catalogue for saw blade range.	KX-600	0 - 13000	
---------------------	---	--------	-----------	--

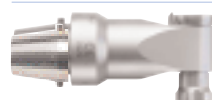

12300	Mini sagittal	NX-600	0 - 25000	
14400 ¹⁾	Mini sagittal	MX-600	0 - 12000	

See Stericut™ catalogue for saw blade range.

¹⁾ MPX only ²⁾ MDX only

See page 13 for other chuck drive options.

Order code	Final drive	Type	Speed rpm
------------	-------------	------	-----------

Minimal invasive acetabular reamer attachment

12100	Accepts De Soutter acetabular reamer shells	MI-600	0 - 150
-------	---	--------	---------

See page13 for range of acetabular reamer shells



Radiolucent drill adaptor

11980	Accepts A.O. radiolucent drill attachment	TX-600	0 - 950
-------	---	--------	---------



Pulse Lavage attachment

11990	Accepts P10 series of pump sets	LX-600	0 - 1700 pulses/min.
-------	---------------------------------	--------	----------------------

See page 11 for range of P10 pump sets



Wire/pin driving attachment

11930	0.6 - 3.2 mm capacity	WX-600	0 - 950
-------	-----------------------	--------	---------

See Stericut™ catalogue for wire and guides range.



MDXMultiDrive™

Modular and procedure specific - battery and



MDX-610 Hi-Torque
*fitted with
A.O. quick release
drill attachment and
sterile battery.*

MDX-610 HANDPIECE SPECIFICATION

Order code	Speed	Cannulation	Voltage*
1286444	0 - 950 rpm	4.1 mm	9.6 Vdc



* Cordless 9.6 Vdc battery or corded direct from power supply module.
See pages 7 and 8 for details of power options.
See pages 3 and 4 for attachment range.
See page 12 for sterilisation containers.

POWER OPTIONS



Small sterilisable
battery



Large sterilisable
battery



Small aseptic battery



Large aseptic battery



Corded system

- Superior handpiece reliability with high performance, hermetically sealed, hybrid motor controller.
- Single handed, instant forward, reverse and safe mode on MDX handpiece.
- 2 speed trigger function on saw handpieces.
- State of the art battery charging technology means batteries are charged to optimum capacity.
- Class leading power to weight ratio.

corded system



CDX-600
fitted with saw blade
and aseptic housing
carrying small
battery.



KDX-600
fitted with saw blade
and aseptic housing
carrying large battery.

HANDPIECE SPECIFICATION

Type	Type code	Order code	Speed	Voltage*
Sagittal saw	KDX-600	1282554	0 - 14000 cycles/min	9.6 Vdc
Reciprocating saw	CDX-600	1282664	0 - 9000 cycles/min	9.6 Vdc



* Cordless 9.6 Vdc battery or corded direct from power supply module. See pages 7 and 8 for details of power options.
See Stericut™ catalogue for range of sagittal and reciprocating saw blades.
See page 12 for sterilisation containers.

ElectricPowerOptions

Ultimate versatility with a choice of three po

STERILE BATTERY

Contains NiCd cells with a maximum capacity of 1200 mAh.

The packs are specifically designed to withstand flash sterilisation techniques.

Extended autoclaving and drying cycles may reduce the capacity and life of the battery pack.

If flash sterilisation is not permissible (please check your market requirements) then the use of aseptic batteries is strongly recommended.



SB-600
9.6 Vdc NiCd
Battery Pack
Order Code 11630




SB-450
10.8 Vdc NiMH
Battery Pack
Order Code 14500

ASEPTIC BATTERY SYSTEMS

The simple aseptic battery technique ensures the packs are maintained and used in their optimal condition.

Battery management is readily supervised by the theatre staff as the battery charger can be safely placed within the OR.

Large Battery Pack
AB-600
9.6 Vdc NiMH battery
Capacity 2100 mAh
Order Code 11520

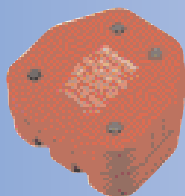

AH-610
Aseptic Housing
Order Code 14430



AS-610
Sterile Shield
Order Code 14420



Small Battery Pack
AB-450
10.8 Vdc NiMH battery
Capacity 1000 mAh
Order Code 14510



AH-450
Aseptic Housing
Order Code 14520



AS-450
Sterile Shield
Order Code 14530




Loading and replacing the aseptic battery pack



Sterile scrub technician places sterile shield onto the aseptic battery housing.



Circulating nurse presses the non sterile battery into the battery housing.

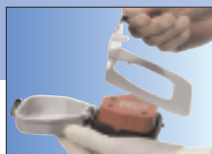

Circulating nurse removes the sterile shield.



Sterile scrub technician closes battery housing lid and secures locking knob.



Sterile scrub technician locks aseptic battery housing into the handpiece.

wer options

BATTERY CHARGER BC-610

Order Code 14750

The BC-610 battery charger allows four batteries to be charged simultaneously.

The unit has a very small footprint and is designed to accept both sterile and aseptic batteries without the use of adaptors.

The charger is easy to operate and provides the user with instantaneous information on a battery's charge status.

A simple display panel allows the user to quickly establish the status of the charging cycle and the general condition of the battery by illuminating a series of indicator bars when it has completed its charging cycle.

Once a battery pack is inserted into the charger it automatically carries out a full diagnostic check and selects the correct charging mode i.e. fast charge, trickle charge etc. The typical re-charging cycle for a used battery is 30 minutes.

Mains power cordsets - 3m long

UK Order Code: 8290

Euro Order Code: 8300

US/Canada Order Code: 8310

Australia/N Zealand Order Code: 8320

Please specify option required.



CA-450
Small battery charging adaptor.
For use with AB-450 battery.
Order Code 14540

CORDED DIRECT POWER

The compact, lightweight corded module simply slides into the bottom of the handpiece and is connected to the power supply via a small diameter flexible cable.

The power supply has two handpiece outlet connections and features an overload indicator which aids the surgeon in achieving optimum cutting performance.

CM-600
Corded module 4m long.
Order Code 11620



PS-600
Power Supply 230 Vac input
Order Code 11610

Power Supply 120 Vac input
Order Code 11660

Mains power cordsets -
3m long


UK Order Code: 8290

Euro Order Code: 8300

US/Canada Order Code: 8310


Australia/N Zealand Order Code: 8320

Please specify option required.

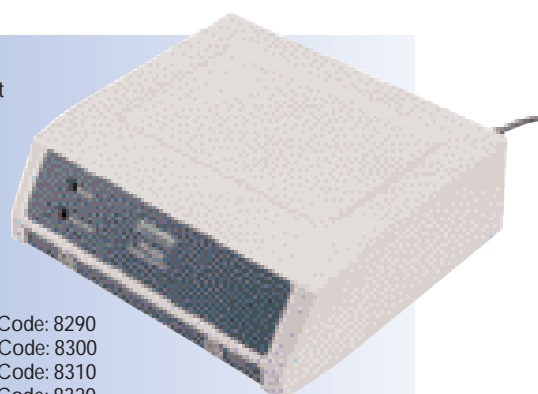

MPX MultiDrive™

Modular and procedure specific - air operate

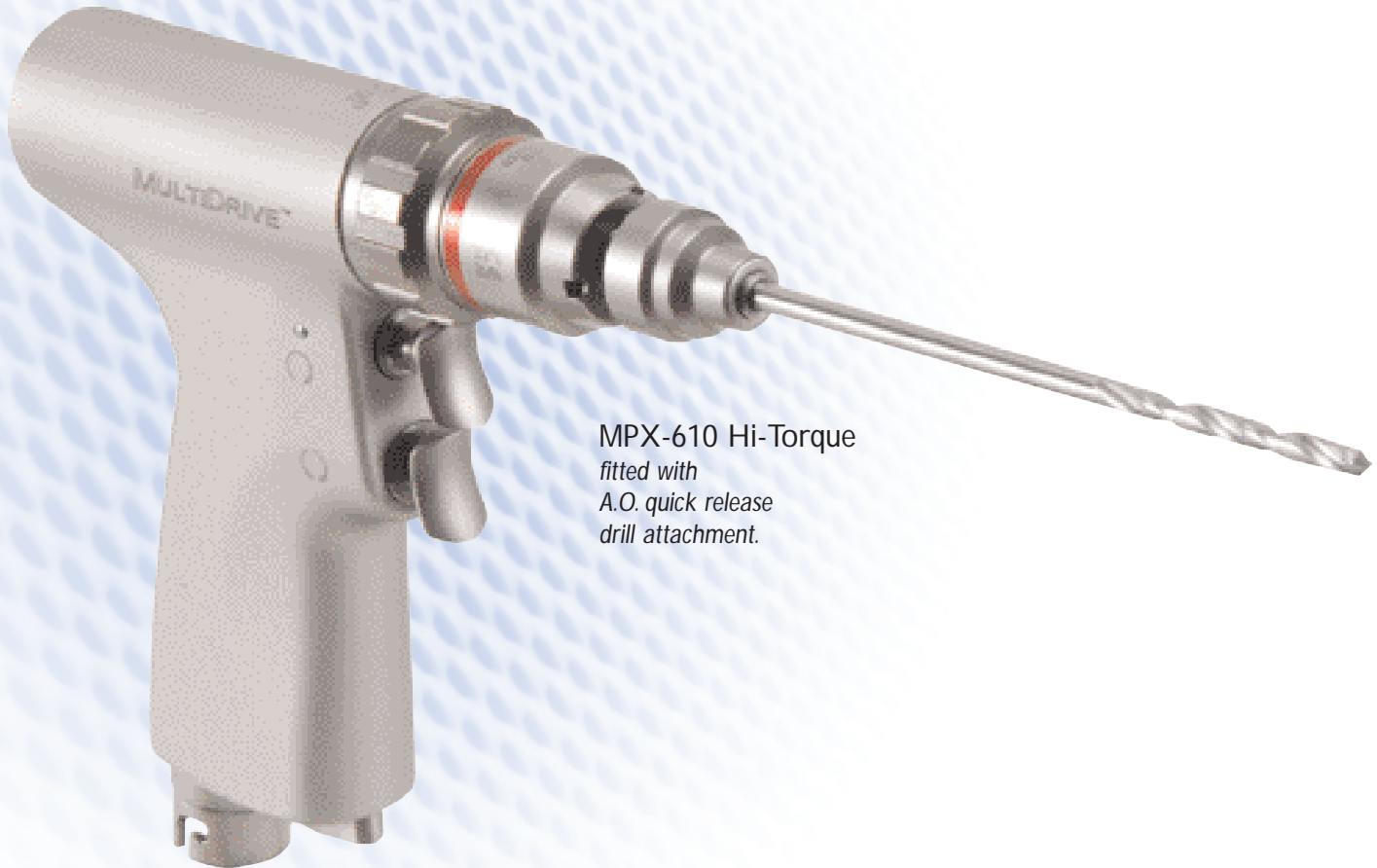

MPX-610 Hi-Torque
fitted with
A.O. quick release
drill attachment.



Three position trigger
function control.
Forward, reverse and safe
operating mode.

MPX-610 HANDPIECE SPECIFICATION

Order code	Speed	Cannulation	Operating pressure
226464	0 - 950 rpm	4.1 mm	6 - 7 bar



ANCILLIARY EQUIPMENT



Modular attachments
See page 3 for details.



Sterilisation container.
See page 12 for details.




Wash Assist hose system.
See page 14 for details.

- Modular versatility for universal trauma and procedure specific handpieces for heavy duty orthopaedics.
- High power, high torque co-axial motor unit
- Single handed, instant forward, reverse.
- Fully machine washable.


d

DPX673
Reciprocating saw

DPX-676
Sagittal saw

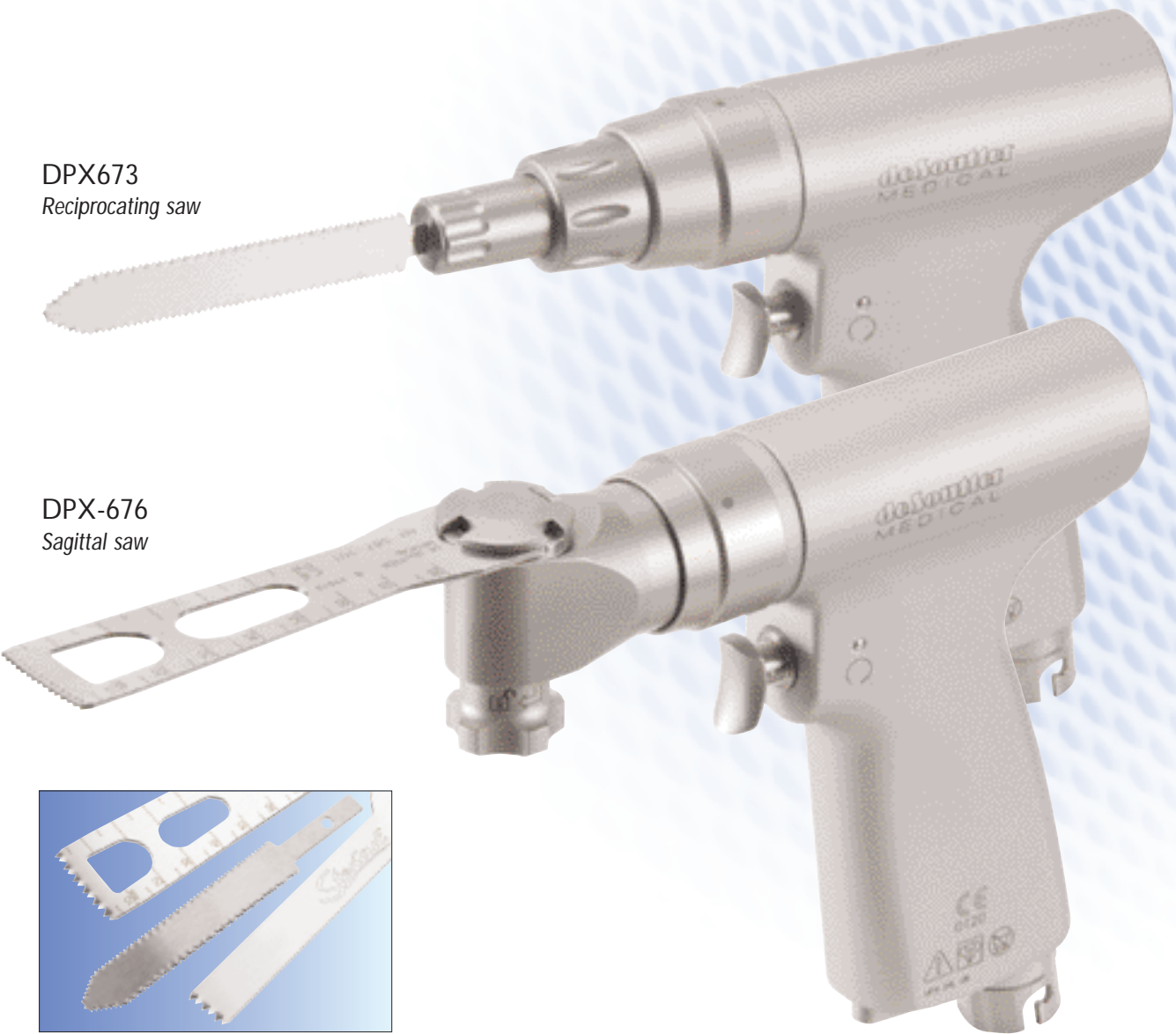


A full range of sterile packed saw blades are available in the De Soutter Medical Stericut™ range. See references below.

HANDPIECE SPECIFICATION					
Type	Type code	Order code	Speed	Operating pressure	
Sagittal saw	DPX-676	226384	0 - 14000 cycles/min	6 - 7 bar	
Reciprocating saw	DPX-673	226394	0 - 11000 cycles/min	6 - 7 bar	

See Stericut™ catalogue for range of sagittal and reciprocating saw blades.
See page 12 for details of sterilisation container.

PLXEcoPulse™

Pulsed irrigation system

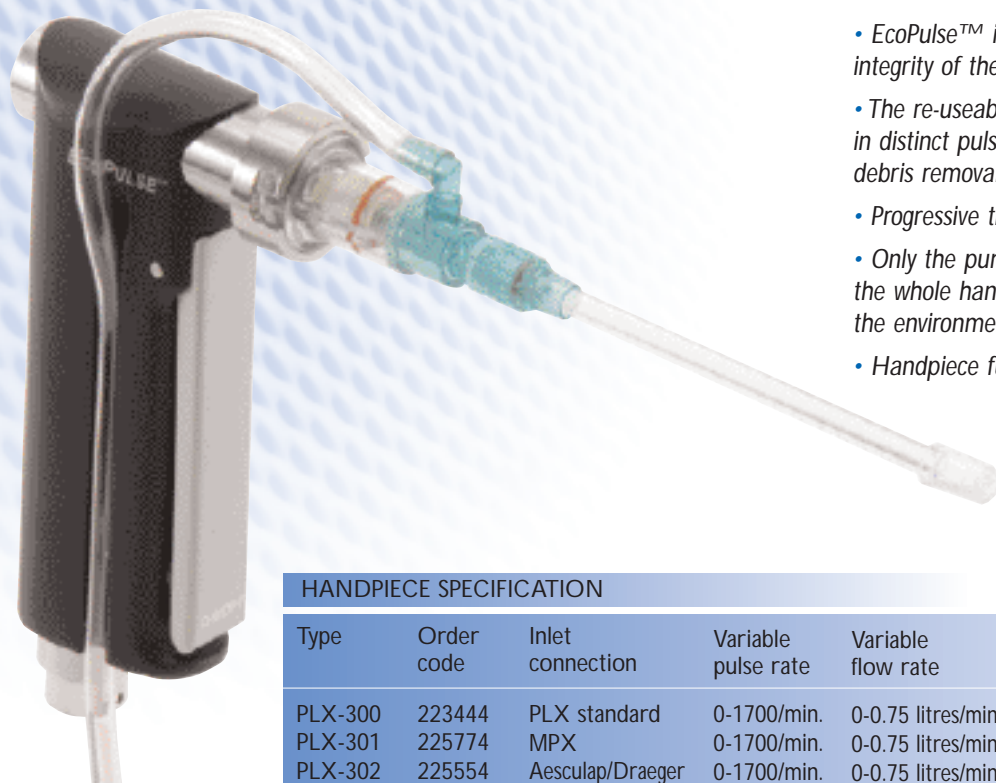

- EcoPulse™ improves the strength and integrity of the bone/cement interface.
- The re-useable handpiece pumps fluid in distinct pulses for effective blood and debris removal.
- Progressive trigger operation.
- Only the pump kits are disposable, not the whole handpiece - saving cost and the environment.
- Handpiece fully autoclavable.




HANDPIECE SPECIFICATION


Type	Order code	Inlet connection	Variable pulse rate	Variable flow rate	Operating pressure
PLX-300	223444	PLX standard	0-1700/min.	0-0.75 litres/min.	6 - 7 bar
PLX-301	225774	MPX	0-1700/min.	0-0.75 litres/min.	6 - 7 bar
PLX-302	225554	Aesculap/Draeger	0-1700/min.	0-0.75 litres/min.	6 - 7 bar
PLX-304	225884	A.O.	0-1700/min.	0-0.75 litres/min.	6 - 7 bar

See page14 for air hose options.


SPLASH SHIELDS

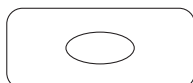

Order code P10-104

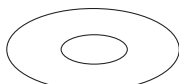


Order code P10-103

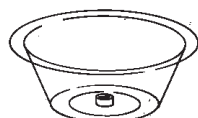

Order code P10-108

All sterile packed pump sets, splash shields and nozzles are supplied 5 per box except P10-108 which is supplied 4 per box.

PUMP SETS

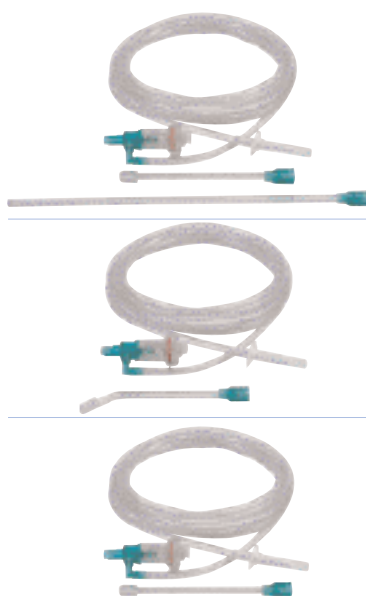


Type	Order code	Nozzle supplied/order code
Hip arthroplasty	P10-101	10 cm shower head P10-105 25 cm femoral shaft P10-107
Knee arthroplasty	P10-102	10 cm angled shower head P10-106
General	P10-100	10 cm shower head P10-105

All pump sets are supplied with pump cartridge, 3 metres of tubing and bag spike connector.

Sterilisation containers

MDX, CDX AND KDX

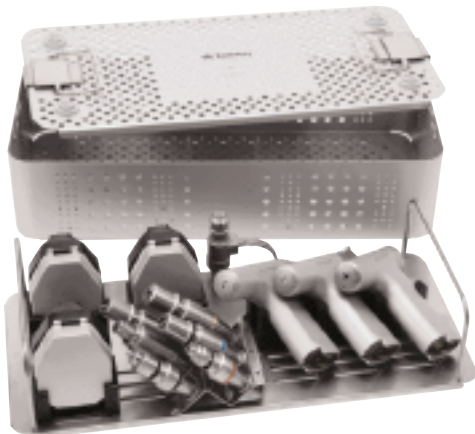

Sterilisation case for system type	Materials	Complete case Order code	Insert only Order code
Aseptic battery 3 handpieces ¹⁾	Aluminium base, lid and stainless steel insert	12170	12180
Aseptic battery 2 handpieces ²⁾	Aluminium base, lid and stainless steel insert	14440	14450
Sterile battery 3 handpieces ³⁾	Aluminium base, lid and stainless steel	11670	11690

¹⁾ Accommodates MDX-610 handpiece, seven attachments, CDX-600, KDX-600, three aseptic battery housings.

²⁾ Accommodates MDX-610 handpiece, seven attachments, KDX-600, two aseptic battery housings and two sterile shields.

³⁾ Accommodates MDX-610 handpiece, seven attachments, CDX-600, KDX-600 and three sterile battery packs.

MPX



Materials	Complete case Order code	Insert only Order code
¹⁾ Aluminium base, aluminium lid and stainless steel insert	14640	14650
²⁾ Aluminium base, aluminium lid and stainless steel insert	14630	14660

¹⁾ Accommodates one MPX handpiece, one DPX handpiece, seven attachments and air hose.

²⁾ Accommodates one MPX handpiece, two DPX handpieces, seven attachments and air hose.

DPX AND PLX



Materials	Complete case Order code	Insert only Order code
Aluminium base, lid and insert	9350	631343


Accommodates air hose and other items in general purpose tray insert.

All sterilisation cases are DIN compatible. Size 538 mm x 254 mm x 142 mm

Accessories

Drive adaptors, reamer shells, drill bits and a

DRIVE ADAPTORS

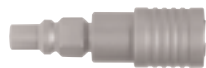

Input drive Hudson Part No.	A.O. Reamer Part No.	A.O. Drill Part No.	Output drive
5440	7970	11550	6.35 mm (1/4") Jacobs chuck
5450	6590	-	Zimmer chuck
5460	6550	-	Harris chuck
-	6560	11580	Hudson chuck
7910	6580	-	Trinkle chuck
5470	7980	-	A.O. Synthes quick release drill chuck
5480	-	-	A.O. Synthes/Protek reamer chuck
5490	6570	-	DIN58809 reamer chuck
11260	7770	-	DHS
-	-	6270	4.0 mm Jacobs Chuck - fits A.O. drill chuck


See Accessory Catalogue publication no. 162 for other drives in the De Soutter Medical range.

ACETABULAR REAMER SHAFT AND SHELLS



A.O. Synthes reamer shaft Part No. 6110

Hudson reamer shaft Part No. 6100

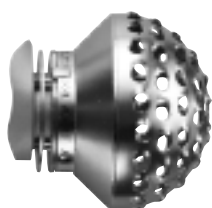

Part No.	Cutting dia mm	Part No.	Cutting dia mm	Part No.	Cutting dia mm	Part No.	Cutting dia mm
5800	40	5880	48	5960	56	6040	64
5810	41	5890	49	5970	57	6050	65
5820	42	5900	50	5980	58	6060	66
5830	43	5910	51	5990	59	6070	67
5840	44	5920	52	6000	60	6080	68
5850	45	5930	53	6010	61	6090	69
5860	46	5940	54	6020	62		
5870	47	5950	55	6030	63		

Fits both standard reamer shaft and minimal invasive reamer attachment.

ir hoses

STERILE PACKED DRILLS

PLAIN SHANK DRILLS




Part No.	Dia. mm	Length mm	Flute length mm
D50-100	1.5	70	22
D50-101	2.0	85	50
D50-102	2.5	95	50
D50-103	2.7	85	50
D50-104	3.2	130	50
D50-105	3.2	180	50
D50-106	3.5	95	50
D50-107	3.5	180	50
D50-108	4.0	156	50
D50-109	4.5	156	50
D50-110	4.5	180	50
D50-111	4.8	127	50
D50-112	5.0	180	50

A.O. DRIVE DRILLS





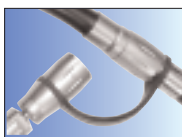







Part No.	Dia. mm	Length mm	Flute length mm	Part No.	Dia. mm	Length mm	Flute length mm
D55-100	1.1	60	20	D55-108	2.7	125	50
D55-101	1.5	85	30	D55-109	3.2	145	50
D55-102	1.5	110	30	D55-110	3.2	195	50
D55-103	2.0	100	50	D55-111	3.5	110	50
D55-104	2.0	125	50	D55-112	3.5	195	50
D55-105	2.5	110	50	D55-113	4.5	145	50
D55-106	2.5	180	50	D55-114	4.5	195	50
D55-107	2.7	100	50	D55-115	6.0	195	50



'WASH-ASSIST' CO-AXIAL AIR HOSES

The Wash Assist co-axial hose system is recommended where machine washing of air hoses is demanded. By releasing the probe end from hose, which is retained on the hose via the rubber strap, the two hose ends can be joined and sealed prior to machine washing.

Probe	Instrument	3 metre hose - Part no.	4 metre hose - Part No.	5 metre hose - Part No.
	MA7 Schrader MPX, DPX PLX	13000	13100	13200
		13600	13700	13800
	A.O. MPX, DPX PLX	13010	13110	13210
		13610	13710	13810
	AGA Inst. MPX, DPX PLX	13020	13120	13220
		13620	13720	13820
	Draeger MPX, DPX PLX	13030	13130	13230
		13630	13730	13830
	Heyer * MPX, DPX PLX	13040	13140	13240
		13640	13740	13840
	Air Liquide * MPX, DPX PLX	13050	13150	13250
		13650	13750	13850
	DIN * MPX, DPX PLX	13060	13160	13260
		13660	13760	13860
	AGA Air * MPX, DPX PLX	13080	13180	13280
		13680	13780	13880
	Flowtech * MPX, DPX PLX	13070	13170	13270
		13670	13770	13870
	NP05 * MPX, DPX PLX	13090	13190	13290
		13690	13790	13890

* Made to special order

Replacement rubber strap Part No. 632393

HANDPIECE WASHING CAPS

Prevents water entering the handpiece during washing prior to autoclaving



Washing cap for MPX and DPX series handpieces.
Part No. 608073



Washing cap for PLX series handpieces.
Part No. 8110



de Soutter
MEDICAL

De Soutter Medical Limited

River Park, Billet Lane, Berkhamsted, Hertfordshire, HP4 1HL. UK
tel: +44 0 1442 860300 fax: +44 0 1442 860333

email: info@de-souttermedical.com web: www.de-soutter.com

publication no. 167-EN-3